## DIRECTIONS FOR USE

### VACUUM MATTRESS EasyFiX PLUS EasyFiX PLUS II

**English:** Directions for use. **Germa EasyFiX PLUS Vacuum Mattress**  
**Deutsch:** Bedienungsanleitung. **Germa EasyFiX PLUS Vakuummatratze**  
**Français:** Mode d'emploi. **Matelas à dépression Germa EasyFiX PLUS**  
**Español:** Instrucciones de uso. **Colchón de vacío EasyFiX PLUS Germa**  
**Svenska:** Bruksanvisning. **Germa EasyFiX PLUS vakuummadrass**  
**Italiano:** Istruzioni per l'uso. **Germa EasyFiX PLUS Materassino A Depressione**  
**Dansk:** Brugsanvisning. **Germa, EasyFiX PLUS vakuummadrass**  
**Nederlands:** Gebruiksaanwijzing. **Germa EasyFiX PLUS Vacuüm Matras**  
**Suomi:** Käyttöohje. **Germa EasyFiX PLUS -tyhjööpatja**  
**Norsk:** Brukerveiledning. **Germa EasyFiX PLUS vakuum-madrass**

**Distributor:**

Ferno Norden Norway AS  
Innlaget 298  
NO-3185 Skoppum  
Norway

Phone: +47 33 03 45 00

Fax: +47 33 03 45 10

E-mail: [fno@fernonorden.com](mailto:fno@fernonorden.com)

[www.fernonorden.com](http://www.fernonorden.com)

**Manufacturer:**

AB Germa  
Industrigatan 54-56  
SE-291 36 Kristianstad  
Sweden

Tel: +46 (0)44 12 30 30

Fax: +46 (0)44 10 31 79

E-mail: [info@germa.se](mailto:info@germa.se)

[www.germa.se](http://www.germa.se)

English:	Directions for use	Page	3
Deutsch:	Bedienungsanleitung	Seite	7
Français:	Mode d'emploi	Page	11
Español:	Instrucciones de uso	Página	15
Svenska:	Bruksanvisning	Sida	19
Italiano:	Istruzioni per l'uso	Pagina	23
Dansk:	Brugsanvisning	Side	27
Nederlands:	Gebruiksaanwijzing	Bladzijde	31
Suomi:	Käyttöohje	Sivu	35
Norsk:	Brukerveiledning	Side	39

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## English Directions for use

### 1. Introduction

The Germa EasyFIX PLUS Mattress is a vacuum mattress. The mattress can also be placed on top of a stretcher/trolley, to be used as a patient mattress.

The Vacuum mattress is principally designed for the immobilisation of patients where bone and spinal injuries are suspected. Deep tissue damage can also be avoided by using vacuum mattress.

The vacuum mattress is fitted with 5 handles on both sides. It is also equipped with five patient restraint straps, each with Fix-Lock, Quick-release plastic buckles.

The Vac mattress should be operated with a minimum of 2 trained operators when carrying the mattress a shorter distance, Additional people may be required for heavier patients and when carrying longer distance. Additional carrying belts can also be used to aid manual handling.

In emergency situation the Vac mattress can be moved or dragged in ergonomically advantageous on the ground using the carry handles or by connecting a carrying belt or rope to the Vac mattress, making sure to avoid rough surfaces or sharp edges on the ground along the route of movement.

The EasyFIX PLUS Mattress consists of an outer shell made from airtight material, which has a number of inner chambers which contain granules. When removing the air from the Vacuum mattress using the supplied vacuum pump, the shell becomes rigid and can be shaped around the patient as the air is removed.

When the air has been fully evacuated from the EasyFIX PLUS Vacuum Mattress, it is rigid and it is not possible to mould it any more. The mattress wraps around the patient offering an optimal support and stabilization, without exerting any pressure on the patient.

The mattress can also be used to carry a patient in water. (See precautions section).

Recommendation: Usable life of 5 years

### 2. Precautions

- Never mould the Vacuum Mattress on the top of the head or under the foot soles as it might cause undesirable pressure against the spine.
- On a regular basis check that the vacuum pump is functioning correctly, and physically test it by deploying the vac mattress by removing air from it until the shell becomes hard.
- Avoid handling the mattress next to sharp objects. If the mattress is being dragged along the ground, be aware of obstacles on the ground, such as sharp stones, sticks, glass, sharp edges etc. Always use intended cover on the valve
- Check pump, Vacuum Mattress and carrying straps for any visible damage or wear and tear. Repair any damage as required before re-using the Vacuum mattress. Change any defect parts immediately. Ensure that all parts are present and operating, including the valve cap, patient restraints and buckles and head strap
- Avoid handling the mattress next to open fire.
- On a current basis check that the mattress has no visible cracks.
- The vacuum mattress is unstable in water. You must never leave a patient unattended in water at any time.
- When transporting in an aircraft, please note that the atmospheric pressure in the cabin may drop when reaching higher altitudes and the vac mattress will therefore become less rigid. In this situation make sure to evacuate more air from the mattress.
- Own modifications on products not allowed – Your medical specialist to be responsible for any modification.
- When decided to take the product out of service and destroy – make sure to follow your local instructions or contact us for additional information.

### 3. Description

The vacuum mattress consists of an outer shell made of hard-wearing PVC coated fabric and a unique granulate-filled inner bag. When the air is evacuated from the vacuum mattress the individual granulates pack up together into a tight mass, making the mattress rigid. Before and while the air is evacuated the mattress may be moulded and positioned according to the shape of the patient, to give best immobilisation. An optimum support of the patient is obtained when all air has been evacuated. This can be determined by the fact that it becomes more difficult to pump as less air remains in the vac mattress, to the point where the pump cannot be moved any more once all the air has been removed. The mattress is equipped with 5 handles on both sides. It is also equipped with five colour coded patient restraint belts with Fix Lock buckles. When not in use, store the belts and fasten them by use of the buckles in the gap between the

row of handles and the mattress.

When the Vac mattress is not being used, and is not required to be used as a mattress on top of a trolley, the vac mattress can be folded and packed to store in the supplied protective bag together with pump and repair kit, taking up less space.

#### 4. Operating the Vacuum Mattress

1. Put the mattress on an even surface near the patient, trying to avoid any sharp objects that may be on the ground. Check that the valve is closed by turning clockwise until fully closed, if needed then twist anti-clockwise ½ a turn to open the valve. Note that the cap should not be removed from the valve at any time during operation, to avoid it getting lost.
2. Position and adjust the mattress on the floor so that there is an even distribution of granules along all areas of the vac mattress, and that it is flat ready to take the patient, with restraint straps unbuckled and not underneath the vac mattress

3. Check that the valve is closed by turning it full clockwise until it feels tight, and then connects the pump.
4. Place the injured patient in the medically desired position on the mattress, utilising standard adopted manual handling procedures, utilising supplementary equipment such as Scoop and slide sheets to aid in safe transfer of the patient.

5. Make sure to support stabilisation of the patient head and neck. Mould the mattress around the patient body to give effective immobilisation, such that it becomes stable and comfortable to the patient. Mould up some of the granular material in the base of the mattress between the legs of the patient. Ensure that the head flaps are positioned effectively around the patient's head and ensure the head strap is attached

6. Connect the patient restraint straps according to the colour code, from legs to chest, and tighten them one by one. X Chest belts are red and black. Cross them over the chest.

Caution! Make sure that the vac mattress is not pressing on top of the patient's head or on the soles of their feet

7. Evacuate the air using the vacuum pump and continue moulding the mattress around the patient's body as more air is removed and as the mattress becomes more rigid. Tighten the restraint straps further as the vac mattress becomes stiffer to give better moulding around the patient. Check that there is no additional pressure to the top of head or under the feet. Check that the unit is stable, give good immobilization and is comfortable for the patient.
8. The patient is now ready to be transported.

#### 4.1 Use of the Vacuum Mattress as a Stretcher/Transfer Mattress

The Vacuum Mattress can be used as a Stretcher/Transfer Mattress, when the air is not removed from the unit.

1. Ensure that the vac mattress valve is opened to allow air to flow into the mattress. Fold the side sections of the mattress over, such that the two side sections meet in the middle and lay flat over the base of the mattress.
2. Preferably place a sheet over the mattress.

3. Place the patient on the mattress and use the patient restraints on the stretcher to secure the patient.

## 5. Cleaning and Disinfecting

After each use the vac mattress should be wiped with a moist cloth and then dried. If dirty the mattress should be washed with soapy water and rinsed with clean water. Allow the mattress to dry completely before storage or usage. If the mattress needs disinfecting first clean the mattress carefully and let it dry. Wipe the complete mattress with a solution containing 70% alcohol or use the disinfectant called LYSETOL AF.

## 6. Maintenance and Repair

At frequent regular intervals the vacuum mattress, valve and pump must be checked for any visible damage.

At the same time the vacuum mattress must be tested for function and the suction operation of the suction pump must be checked according to Section 7.

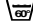
Leaks should be checked for by visual inspection or by slightly filling of the mattress with air, covering it with soapy water all over the surfaces and the welds, then watch for air bubbles that may appear, which will indicate a leak. When filling the mattress with air to check for leaks, ensure that the mattress is not filled so much such that it creates a pressure within the mattress, which may then cause the mattress edge weld to split.

If any leak is found on the mattress, first clean the area around the leak and allow it to dry.

Find a repair patch of a suitable size from the repair kit and then spread the supplied glue all over the patch and all around the area of the vacuum mattress where the patch will be applied.

Wait for the glue to become "dust dry" (approx. 2-3 minutes) and then place the patch on the area of glue and leakage. Apply pressure to the whole area of the patch; making sure good contact has been established in the whole area and that any air bubbles are removed.

The repair will have full strength after 24 hours, but can be used within a few minutes in case of emergency

Note: Patient restraint straps are detachable, washable  and replaceable, washing in laundry bag, see label at the foot end restraint. ! If blood gets on the straps, soak in cold water first before washing

## 7. Testing the Vacuum Mattress and Pump

Straighten out the mattress on the floor, close the valve by turning it fully clockwise and connect the pump. Evacuate the air until the mattress becomes rigid and leave it for about 2 hours. Check to ensure that the mattress is still as rigid. Any loss of rigidity may indicate a leak in the vacuum mattress or through the valve. Check both for leaks and correct functioning of the valve cap.

See Section 6 Maintenance and Repair or check the valve for leakage from the vac mattress. If OK, try to insert a new valve and repeat the test again.

If the pump does not seem to be removing air from the vac mattress, check the pump function by removing it from the vacuum mattress, place a finger over the orifice at the end of the pump tubing and pump a few times. A negative pressure should be created in the tubing and should be felt at your fingertip.

If the pump does not function do not try to repair it, return it for service to your local dealer or to a Germa service centre.

## 8. Storage

After each use, check to ensure that the mattress is still operating correctly and made ready for its next use. Open the valve by turning it fully clockwise to fully close then ½-1 turn anti-clockwise to open and put the mattress onto an even surface. Make sure the granules have been distributed evenly, close the valve by turning clockwise and then fold the mattress in its initial position.

The vacuum mattress can be stored on top of the stretcher, or folded in the storage bag, such that it will fit in the designated storage area.

## 9. Specifications

The Germa EasyFix PLUS Vacuum Mattress is tested in accordance with EN 1865 and EN 1789  
The Germa EasyFix PLUS Vacuum Mattress is in conformity with Council MDR 2017/745 Class 1 concerning Medical Devices

### Dimensions:

Length x Width x Height: see product information below  
Weight: see product information below  
Heat resistance: 70°C  
Cold resistance: -30°C

Capacity for foot pump: The pump has capacity to empty the mattress according to EN 1865

### Materials:

Part	Material
Top and bottom material:	Polyester fibres PVC coated (flame resistant)
Extra reinforcement material:	Polyester fibres PVC coated (flame resistant)
Innerbag:	Polyester
Wood:	Birch

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**Report no. :**

(To be assigned in by  
Manufacturer)

**Please print or type or fill in by MS Word, all sections a. – i. Fax or e-mail to your supplier. Please use repair order when reporting repairs not associated with complaints.**

For claims please specify :  
Purchase Date

## Description of product

**Description of Incident/Problem:** 1) What, 2) How?

When? In connection with: 1) transport, 2) unpacking, 3) preparing, 4) test, 5) cleaning, 6) normal use, 7) training, 8) disassembly, 9) assembly, 10) storage, 11) other:

**Info on the product, prior to problem,**

Storage: 1) stock, 2) site of use, 3) packed, 4) unpacked, 5) softpack, 6) other:

Used cleaning and disinfecting agents/Sterilizing method:

Pattern of use: 1) stand-by for emergency, 2) routine use, 3) stand by, 4) other:

Environment: 1) Indoor, 2) Outdoor, 3) Heat, 4) Cold, 5) Humidity, 6) Sunlight, 7) Gases, 8) other:

Which **actions** have been: a) taken, b) planned by the reporter or user at the site of incident?

1) examination of the product on site, 2) replacement of the product, 3) other:

---

Suggested actions by the distributor/complainant

---

---

**Returning of samples/faulty product(s),**

- 1) Already shipped/Will be shipped to \_\_\_\_\_,
- 2) Pending request from manufacturer,
- 3) Product/parts not available:

---

Date:

---

**We acknowledge the receipt of this report by manufacturer and will start processing it immediately. Please use the assigned complaint report number for further reference.**

**Your contact person in this matter is:** \_\_\_\_\_

Name:

Direct Tel.

Direct Fax.

Init/Date:

---