



# FERNO® XTRIC8 RESCUE SLED



The Ferno XTRIC8 Rescue Sled is a versatile and cost-effective solution for universal rescues including lowland, mountain, confined space and high angle extrications.

The XTRIC8 is a compact and durable system. Smart features decrease the time and effort necessary to package patients while providing traditional functionality and tactical evacuation methods in an easy to use device.

## Key Features

- XTRIC8 is supplied as a complete, ready to use unit with all necessary hardware
- The XTRIC8 is pre-rigged for immediate deployment out of the bag
- Includes an internal patient harness for maximum patient security during vertical lifts
- Colour coded, 4-point lifting bridle is pre-set to aid faster implementation
- Reflective external webbing for enhanced visibility in low light conditions
- Straps and buckles are glove-friendly and designed for efficient attachment and release
- Robust carry bag with reflective webbing, four large pockets and adjustable backpack harness
- XTRIC8 is intuitive to use and minimises user training

### XTRIC8 Load Rating



350 KG



250 KG

Internal patient harness keeps patients secure and prevents unwanted movement during vertical rescues



Glove-friendly buckles and quick release tabs improve operational efficiency



XTRIC8 lifting bridle is colour-coded and labelled to enable faster deployment




Reinforced lifting points are colour-coded to match the XTRIC8 lifting bridle for easy identification in the field



Six large rubber handles are glove-friendly and evenly spaced around the XTRIC8 for even load distribution

Drag handle and all XTRIC8 external webbing is reflective for low light visibility

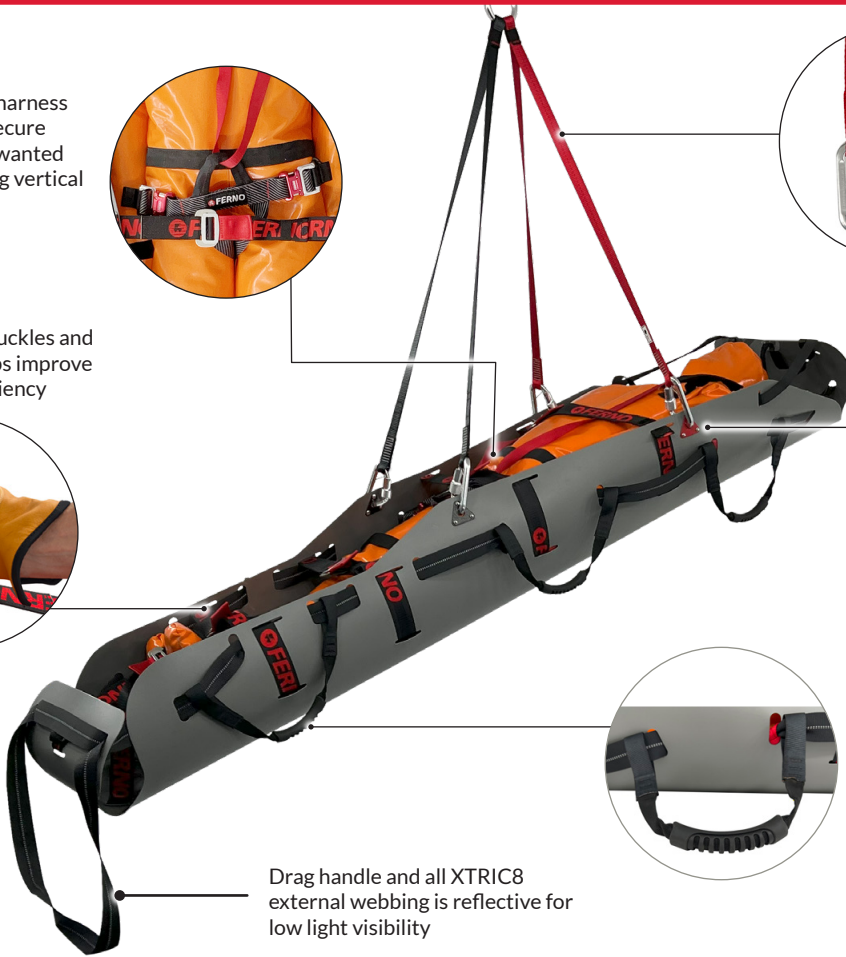

### XTRIC8 RESCUE SLED

<b>Product Code</b>	FWEXTRIC8
<b>Inclusions</b>	XTRIC8 Rescue Sled pre-rigged with internal patient harness assembly, external restraints, external straps and handles
	1 x Maillon Rapide for vertical lifting
	1 x Colour coded 4-Point Lifting Bridle
	1 x Storage Bag with backpack straps
<b>Length</b>	2540 mm
<b>Width - Widest Point</b>	914 mm
<b>Weight - XTRIC8 only (rigged)</b>	7 kg
<b>Weight - XTRIC8 (rigged) and Bag</b>	8.5 kg
<b>Load Limit - XTRIC8</b>	350 kg - HORIZONTAL / 250 kg - VERTICAL

### XTRIC8 LIFTING BRIDLE

<b>Weight - Lifting Bridle</b>	0.9 kg
<b>Load Limit - Lifting Bridle</b>	500 kg
<b>Bridle Karabiners</b>	VYO N-2853G - Aluminium Alloy (SAE 7075) Triple Action Karabiner
<b>Bridle Karabiner Dimensions</b>	124 mm x 80 mm / Gate Opening: 22 mm
<b>Bridle Karabiner Ratings</b>	Minimum Tensile Strength: 23 kN / Minimum Gate Strength: 16 kN
<b>Bridle Karabiner Standards</b>	ANSI Z359.12-2019; CSA Z259.12-16; EN 362:2004

- Ferno and manufacturers of components reserve the right to change specifications without notice.
- General specifications are rounded to the nearest whole number. Metric conversions are calculated before rounding the English measurements. For more information contact Ferno Customer Service or your local distributor.